


Tomato ^{Super} CAST™

Orthopedic Casting Tape

CAST / SPLINT
2in1 SAFE PAD / COTTON PAD
ELASTIC BANDAGE
STOCKINET




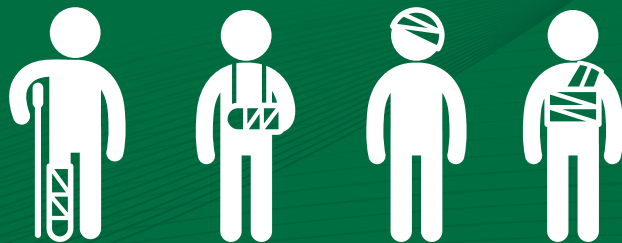


TOP MANAGEMENT FOR TOMORROW of Fiber Cast and Splint

Tomato M&C Co., Ltd. was established in February 2005 and became a expert of manufacturer that is specializing in synthetic orthopedic cast and splint.

We, Tomato M&C are trying to establish the identification for tomato casts and tomato splint through friendly product designs and steady upgrades of product core materials. We are distributing our qualified products and services to **domestic and overseas customers around 30 countries such as the United States, Japan, China, Southeast Asia, MENA, Europe, South America.** Tomato M&C Co., Ltd. as a one of leading manufacturers, we are supplying **Tomato brand products, OEM products, production facilities, and raw materials** under **ISO13485, ISO14001, ISO9001, USA-FDA, CE DOC, KGMP** quality management system. We are the only one manufacturer that produce products with a one-stop button operating system through automatic packaging facilities for synthetic resin casts in Korea and neighboring countries.

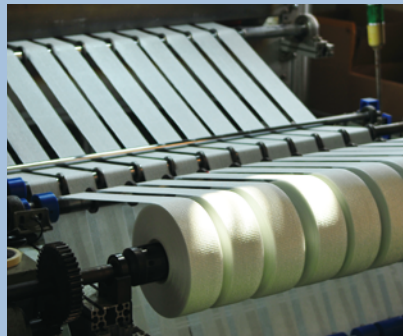

Conceptual diagram of tomato cast production facilities (for Tomato Fiber Cast / Tomato Poly Cast)

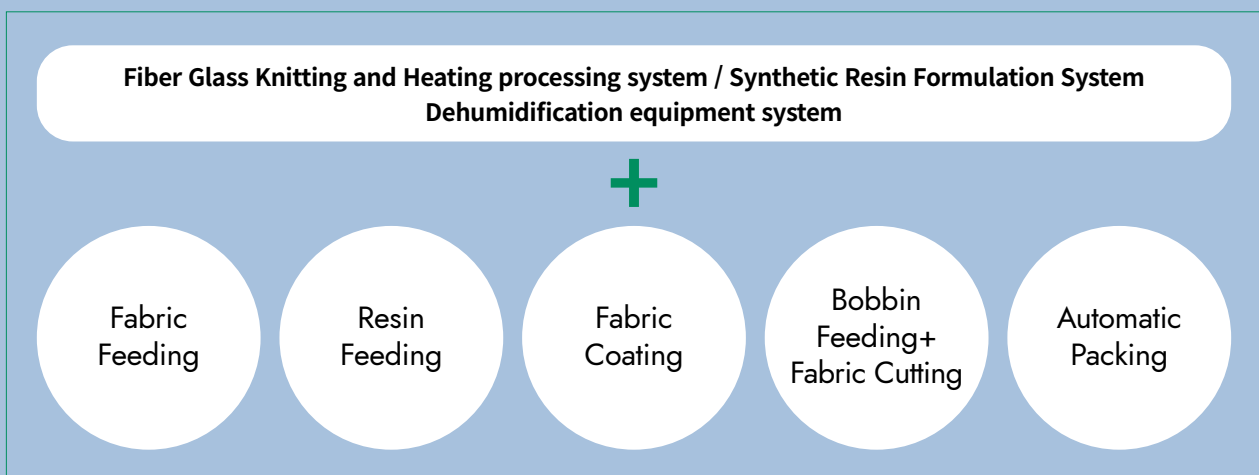

Tomato Cast / Tomato Splint's manufacturing system is one of Tomato M&C's proud sales items, and we have grown with the Korean synthetic orthopedic fiber cast and splint industry. Now we have the most efficient manufacturing system for production and maintenance.

Tomato ^{Super} CAST™

FEATURES

- Comfortability: Smooth unrolling offers doctors can wrap it with less pressure and patient will have no pains during the applying procedure.
- Conformability and mouldability : Our own produced extreme soft fabric makes it easier to wrap and shape the contours of the human body.
- Durability for safety: Tomato Cast is strong to impact and durable on the edges.
- Tack-free and no layer separation
- Excellent X-ray translucency

APPLICATION INSTRUCTIONS


Please refer to our websites (www.tomatomnc.com) for instruction video and more details.




- 1 Slip the two layers of Tomato Safe Pad on the affected body part.



- 2 Gloves must be worn.

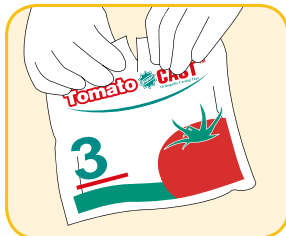


- 3 Open only one roll of Tomato Cast at a time.

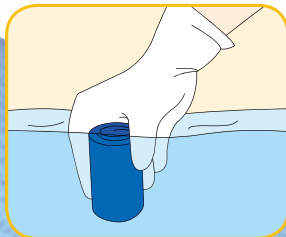

- 4 Immerse the roll in the water in a temperature of 20~24 degrees for 3~5 seconds and squeeze 3~4 times in the water for the complete water penetration into the body of the roll. This process will help strengthen adhesion and make it easy to come off the roll.
* Do not put it in the water for too long.



- 5 Wrap spirally, taking care to avoid excessive tightness.
* Finish wrapping it up within 1 minute and 10 seconds after the pouch is open.

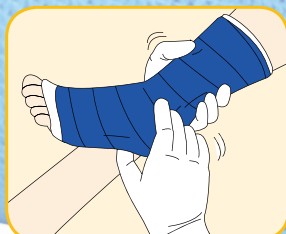

- 6 Smooth and rub the surface to achieve good contact between layers. In 3-5 minutes after immersion in the water, Tomato Cast is sufficiently cured for further molding.

PRODUCT INFORMATION

MODEL NO.	SIZE
TRC-001F/TRC-001P	1in x 2yds (2.5cm x 1.8m)
TRC-002F/TRC-002P	2in x 4yds (5.0cm x 3.6m)
TRC-003F/TRC-003P	3in x 4yds (7.5cm x 3.6m)
TRC-004F/TRC-004P	4in x 4yds (10.0cm x 3.6m)
TRC-005F/TRC-005P	5in x 4yds (12.5cm x 3.6m)
TRC-006F/TRC-006P	6in x 4yds (15.0cm x 3.6m)

COLOR INDEX



※ The color shown in the picture may be different from the actual color.

ORDERING CODE

TRC - 00XX - XX

COLOR INDEX
F: FIBERGLASS / P: POLYESTER
WIDTH
TRC: TOMATO RIGID CAST TSC: TOMATO SOFT CAST

* e.g. - TRC - 003F - WH: tomato rigid cast, 3 inch, fiberglass, white

Removal: Tomato Cast can be removed with a traditional cast saw. Care should be taken to avoid injury to the patient when a small amount of underpadding has been used.

INTENTION TO USE

Tomato Cast has been designed for immobilization of fractures, sprains and strains, ect.

PRECAUTION

1. Each roll should be gently squeezed prior to opening to check suitability for application. The unopened roll should be soft and flexible, Do not use the roll that feels hard.
2. The polyurethane resin could adhere firmly to unprotected skin, hair and clothing, which may cause a mild irritation. Thus, always wear protective gloves to apply and use alcohol or acetone before to remove the polymer.
3. An uncomfortable increase of heat can be caused during hardening by excessive layering of tape or immersing in water at temperature above 27°C.
4. All application should be followed instructions and directions of a doctor or qualified medical professional.
5. Products should be stored out of reach of children.



Tomato SOFT CAST™

FEATURES

- Easy to be cut off by bandage scissors
- Enough fixation for support
- Skin tight applications without padding
- Shoes can be worn with the cast on
- Suitable for a child in cast removal



INTENTION TO USE

Tomato Soft Cast is intended for secondary casting and a comprehensive wrap to control swelling. Tomato Soft Cast, more flexible casting material is typically used for functional immobilizations that can tolerate some movement (ex. minor broken bones, or more casual strains, sprains, and **ligament damage**).

PRODUCT CHARACTERISTICS

Tomato Soft Cast consists of a knitted polyester fabric impregnated with a polyurethane resin. Exposure to moisture or water initiates a chemical reaction which causes the roll to set. The finished wrap will retain its shape while remaining elastic to the touch

PRODUCT INFORMATION

MODEL NO.	DESCRIPTION	SIZE
TSC-002P	Tomato Soft Polyester Cast	2in × 4yds (5cm × 3.6m)
TSC-003P	Tomato Soft Polyester Cast	3in × 4yds (7.5cm × 3.6m)
TSC-004P	Tomato Soft Polyester Cast	4in × 4yds (10cm × 3.6m)
TSC-005P	Tomato Soft Polyester Cast	5in × 4yds (12.5cm × 3.6m)


Tomato ^{Super} SPLINT™

FEATURES

- Comfortability & hygiene : Instead of soaking in the water, as water spray is possible on the mesh part, it gives more pleasant feeling to the patient because there is no water on patient side.
- Conformability & mouldability: Tomato Splint is very easy to shape the contours of human body.


APPLICATION INSTRUCTIONS

Please refer to our website (www.tomatomnc.com) for instruction video and more details.

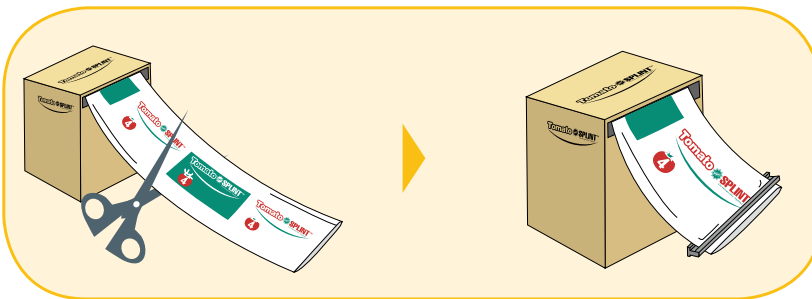

1A Take out Tomato Splint covered with long sealed foil pouch from the box, measure the length to apply and cut off with a bandage scissors.


1A-1 Seal up immediately with a clip (sealing bar) provided inside each box to prevent air or moisture from coming inside the package.

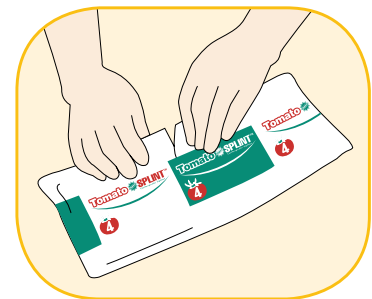

2 Choose the desired size for the affected body part. Generally no gloves are necessary.




3 Except using traditional methods, you can spray room temperature water on the out side part (Green mesh part) of Tomato Splint.

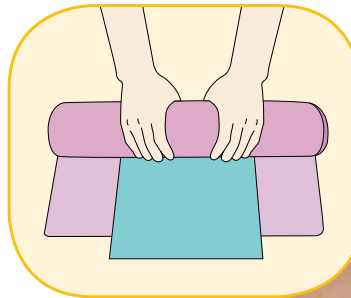

4 Get ride of excessive water on the splint after rolling up with the dry towel. Depending on times, it is not needed to get ride of water totally. In case of need, just tap on the splint surface with dry towel.



5 Wrap Tomato Splint with a roll of elastic bandage to secure and mold it freely.

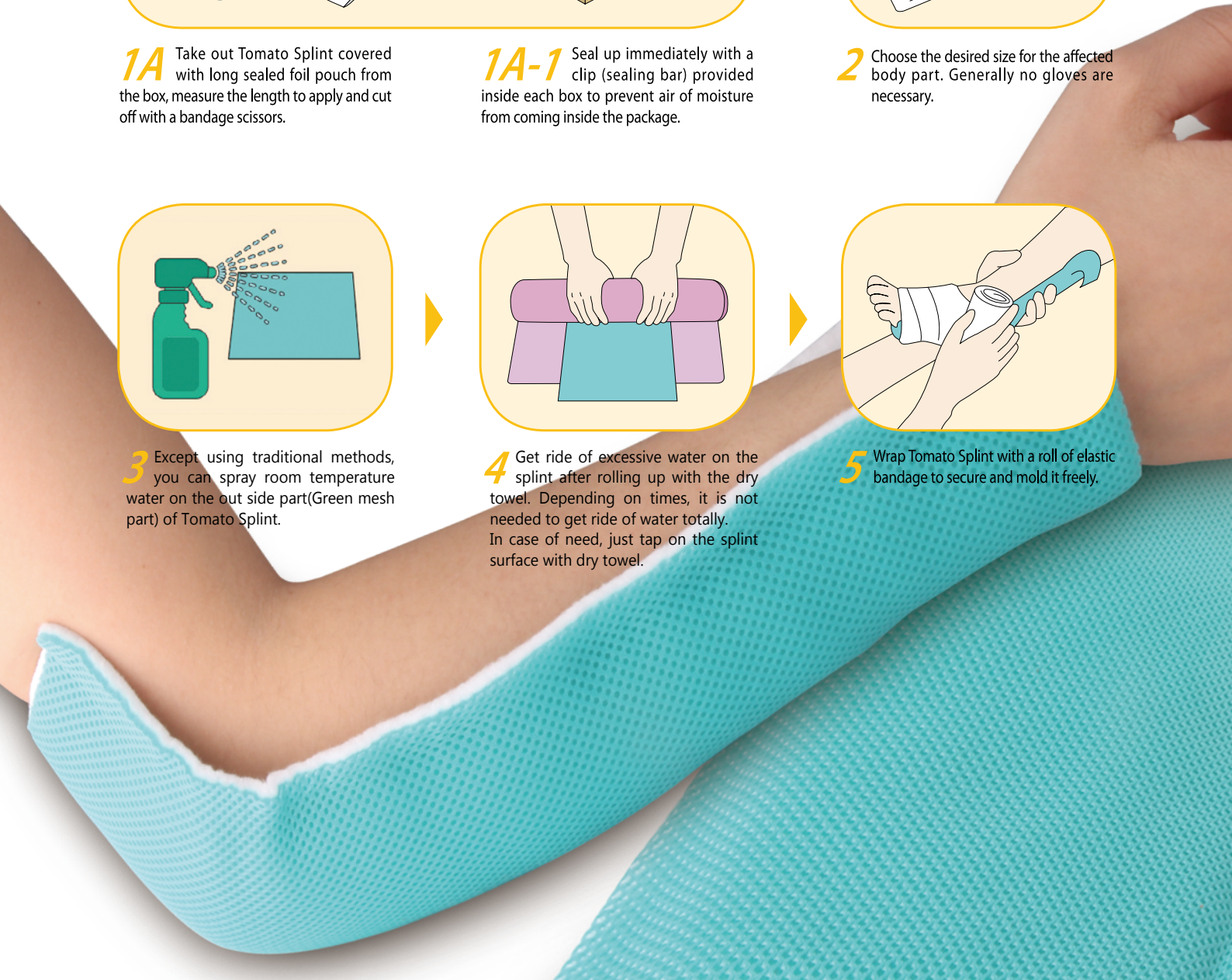


INTENTION TO USE

Tomato Splint has been designed for use in emergency treatment for fractures, sprains and strains. Common practices are, after fractures, swelling would cause a cast to function improperly. Tomato Splint can be used until swelling has been decreased and ready to apply complete cast.

PRODUCT INFORMATION

MODEL NO.	TYPE	SIZE
TPS-210F / TPS-210N	PRE-CUT	2in×10in(5.0cm×25.0cm)
TPS-312F / TPS-312N	PRE-CUT	3in×12in(7.5cm×30.0cm)
TPS-335F / TPS-335N	PRE-CUT	3in×35in(7.5cm×87.5cm)
TPS-415F / TPS-415N	PRE-CUT	4in×15in(10.0cm×37.5cm)
TPS-430F / TPS-430N	PRE-CUT	4in×30in(10.0cm×75.0cm)
TPS-530F / TPS-530N	PRE-CUT	5in×30in(12.5cm×75.0cm)
TPS-545F / TPS-545N	PRE-CUT	5in×45in(12.5cm×112.5cm)
TPS-630F / TPS-630N	PRE-CUT	6in×30in(15.0cm×75.0cm)
TPS-645F / TPS-645N	PRE-CUT	6in×45in(15.0cm×112.5cm)
TRS-002F / TRS-002N	ROLL	2in×5yds(5.0cm×450cm)
TRS-003F / TRS-003N	ROLL	3in×5yds(7.5cm×450cm)
TRS-004F / TRS-004N	ROLL	4in×5yds(10.0cm×450cm)
TRS-005F / TRS-005N	ROLL	5in×5yds(12.5cm×450cm)
TRS-006F / TRS-006N	ROLL	6in×5yds(15.0cm×450cm)
TRS-008F / TRS-008N	ROLL	8in x 5yds (20.0cm x 450cm)

ORDERING CODE

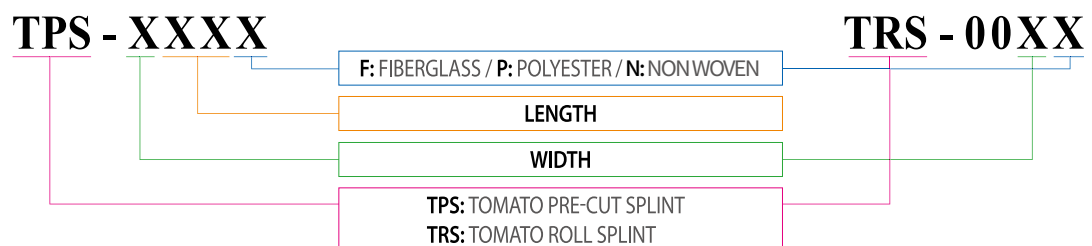


* e.g. - TPS - 312P: Tomato pre-cut type splint, 3 in × 12 in, polyester

* e.g. - TRS - 004F: Tomato roll type splint, 4 inch, fiberglass

EVAPORATION DURING THE HARDENING PROCESS

This picture evidences that our splint evaporates during the hardening process. This makes patients' condition most comfortable.

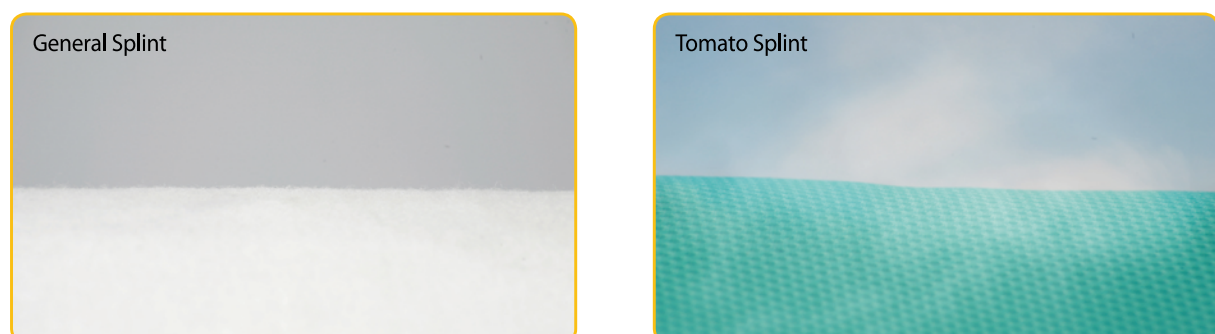

※ Test Condition : Temperature -10±1℃

Tomato 2in1 SAFE PAD™

FEATURES

- **Comfortability**

Tomato  2in1 SAFE PAD™ reduces the water absorption and dries quickly to keep a patient more comfortable.

- **Excellent Air Permeability**

Tomato  2in1 SAFE PAD™ is well ventilated and reduces skin irritation caused by sweat and moisture.

- **Superior Cushioning**

Tomato  2in1 SAFE PAD™ cushions the impact effectively.

PRODUCT CHARACTERISTICS

- Tomato  2in1 SAFE PAD™ is the product mixed with PAD and STOCKINET and consists of polyester fabric.
- Tomato  2in1 SAFE PAD™ is easy to use and reduces the medical treatment time.
- Tomato  2in1 SAFE PAD™ protects the affected area and skin when CAST is removed by CAST SAW.

PRECAUTION

- Tomato  2in1 SAFE PAD™ reduces the water absorption, but avoid taking a sauna or shower without a doctor's permission.
- Excessive actions can lead to a worse injury.

PRODUCT INFORMATION



MODEL NO.	SIZE
TRP-002	2in × 10m
TRP-003	3in × 10m
TRP-004	4in × 10m
TRP-005	5in × 10m
TRP-006	6in × 10m

Tomato PAD™

FEATURES

- Extreme softness with volume makes superior comfortability
- Once binding up, won't be come off easily
- 100% Cotton minimizes skin troubles

INTENTION TO USE

Tomato  PAD™ is intended as a soft sores-preventing inter-layer under solid casting tape.



MODEL NO.	SIZE
TMP-002	Orthopedic pad 2"×2.7m, 100% cotton
TMP-003	Orthopedic pad 3"×2.7m, 100% cotton
TMP-004	Orthopedic pad 4"×2.7m, 100% cotton
TMP-006	Orthopedic pad 6"×2.7m, 100% cotton
TMP-008	Orthopedic pad 8"×2.7m, 100% cotton

Tomato STOCKINET™

FEATURES

- Soft & moderate elasticity offers more comfortability to patient
- It minimizes skin trouble because it is made of cotton

INTENTION TO USE

Tomato STOCKINET is intended as a soft sores-preventing inter-layer under solid cast materials.



MODEL NO.	SIZE
TMS-002	100% cotton stockinette 2"×10m
TMS-003	100% cotton stockinette 3"×10m
TMS-004	100% cotton stockinette 4"×10m
TMS-006	100% cotton stockinette 6"×10m
TMS-008	100% cotton stockinette 8"×10m
TMS-010	100% cotton stockinette 10"×10m

Tomato ELASTIC BANDAGE™

FEATURES

- You can choose either 100% Cotton or 95% Cotton and 5% Rubber case by case.
- Providing breathing comfort in case of bandaging of torso thanks to its elasticity
- Allowing unconstrained movement after bandaging of knee or elbow joint.
- Good value for the money as the amount of usage is less

INTENTION TO USE

Tomato ELASTIC BANDAGE is intended as holding the splint in place and to protect it.



MODEL NO.	SIZE
TME-002R / TME-002C	5cm×2.15m
TME-003R / TME-003C	7.5cm×2.15m
TME-004R / TME-004C	10cm×2.15m
TME-006R / TME-006C	15cm×2.15m