

Tomato  **CAST**™
Orthopedic Casting Tape

CAST
CAST ROLL
STOCKINET
SPLINT



★ US FDA APPROVED ✓

★ 30+ COUNTRIES

★ MADE IN KOREA




TOP MANAGEMENT FOR TOMORROW of Fiber Cast and Splint

Tomato M&C Co., Ltd. was established in February 2005 and became a expert of manufacturer that is specializing in synthetic orthopedic cast and splint.

We Tomato M&C are an established company for Tomato Casts and Tomato Splint through friendly product designs and steady upgrades of product core materials. We are distributing our qualified products and services to **domestic and overseas customers around 30 countries such as the United States, Japan, China, Southeast Asia, MENA, Europe, South America.** Tomato M&C Co., Ltd. as a one of leading manufacturers, we are supplying **Tomato brand products, OEM products, production facilities, and raw materials** under **ISO13485, ISO14001, ISO9001, USA-FDA, CE DOC, KGMP** quality management system. We are the only one manufacturer that produce products with a one-stop button operating system through automatic packaging facilities for synthetic resin casts in Korea and neighboring countries.



Conceptual diagram of tomato cast production facilities (for Tomato Fiber Cast / Tomato Poly Cast)

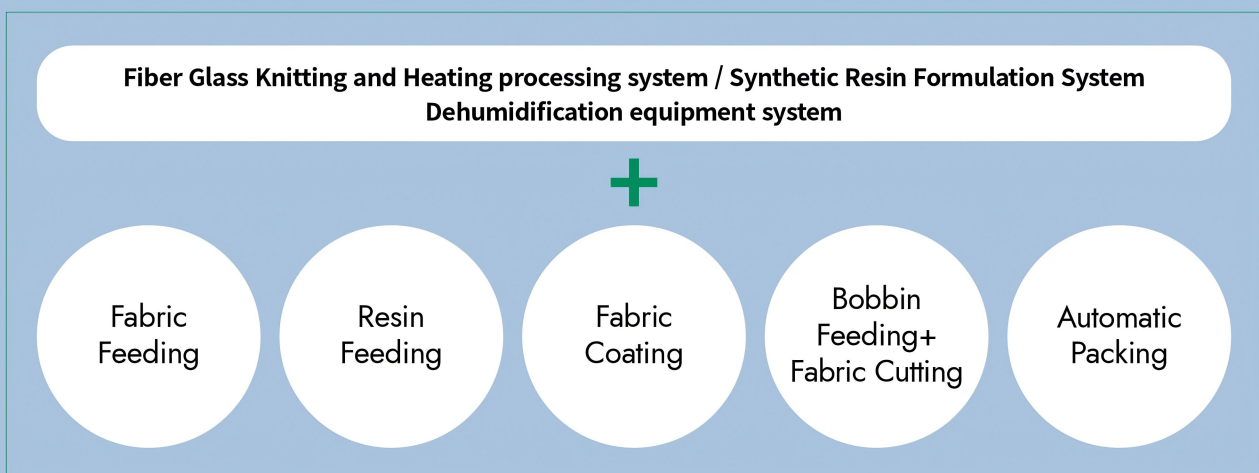


Tomato Cast / Tomato Splint's manufacturing system is one of Tomato M&C's proud sales items, and we have grown with the Korean synthetic orthopedic fiber cast and splint industry. Now we have the most efficient manufacturing system for production and maintenance.

Tomato ^{Super} CAST™

FEATURES

- Comfortability — smooth unrolling lets clinicians wrap with less pressure and no pain to the patient
- Better conformability — molds smoothly around complex joints like ankles and wrists without wrinkling
- Durability for safety — strong on impact and durable on the edges
- Tack-free, smoother lamination — doesn't stick to gloves during application while the layers bond perfectly to each other
- Cleaner edges — no fraying or sharp fiberglass needles to poke and irritate the skin after curing Excellent X-ray translucency

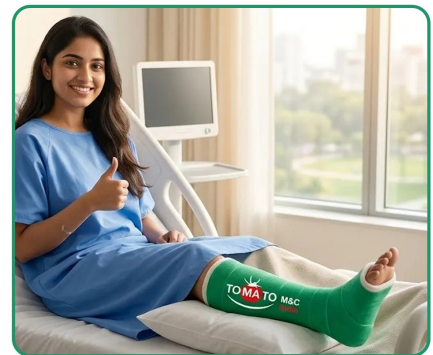

COLOR INDEX



White



Green



Purple



Pink



Blue



Red



Orange



Yellow



Grey



Black



Sky Blue (Blue Clear)



Neon Green



Green Pastel



Pink Pastel




Ocean Green

PRODUCT INFORMATION

REF CODE	WIDTH	LENGTH	TYPE
TRC-002	2 in (5.0 cm)	3.6 m	Fiberglass
TRC-003	3 in (7.5 cm)	3.6 m	Fiberglass
TRC-004	4 in (10.0 cm)	3.6 m	Fiberglass
TRC-005	5 in (12.5 cm)	3.6 m	Fiberglass
TRC-006	6 in (15.0 cm)	3.6 m	Fiberglass


APPLICATION INSTRUCTIONS

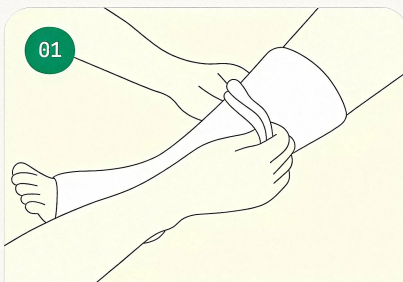

Preparation

STEP / 01

Slip the two layers of Tomato Safe Pad onto the affected body part.

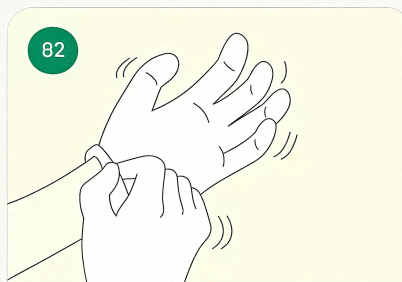

Wear gloves

STEP / 02


Gloves must be worn before handling the cast roll.



Open one roll

STEP / 03

Open only one roll of Tomato Cast at a time to control set time.

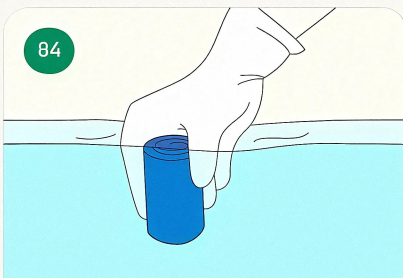

Immerse

STEP / 04


Immerse the roll in water at 20–24°C for 3–5 seconds and squeeze 3–4 times for full penetration. Do not soak too long.



Wrap

STEP / 05

Wrap spirally, avoiding excessive tightness. Finish within 1 minute 10 seconds of opening the pouch.

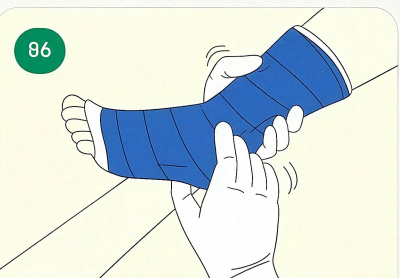

Mould

STEP / 06

Smooth and rub the surface for good contact between layers. Cures for moulding in 3–5 minutes.




FEATURES

- Plush Protective Loft — maximizes shock absorption and cushions friction points to eliminate pressure sores and chafing over bony prominences.
- Advanced Micro-Climate Control — highly breathable fibers rapidly absorb moisture and perspiration, keeping the skin exceptionally dry, calm, and itch-free.
- Precision Conformability — tears cleanly without scissors and self-adheres slightly for a smooth, shift-free wrap around difficult joint contours.
- Optimized Thermal Shielding — provides high-density, uniform protection against the exothermic heat generated during the plaster or fiberglass setting process.
- 100% Premium High-Grade Absorbent Cotton — a soft, fluffy, non-woven fibrous web with uniform density.

PRODUCT INFORMATION

REF CODE	WIDTH	LENGTH	TYPE
SCR-002	2 in (5.0 cm)	3 m	Premium Cotton Blend
SCR-003	3 in (7.5 cm)	3 m	Premium Cotton Blend
SCR-004	4 in (10.0 cm)	3 m	Premium Cotton Blend
SCR-006	6 in (15.0 cm)	3 m	Premium Cotton Blend

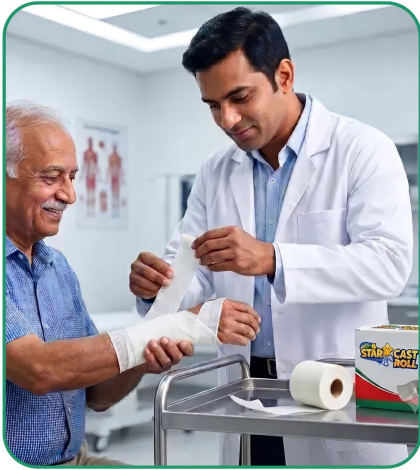



FEATURES

- Conforms smoothly, retains elastic memory
- Resists running, holds shape at cut lines
- Soft, breathable, wicks or absorbs moisture
- Soft & moderate elasticity offers more comfort to the patient
- 100% cotton minimizes skin trouble

PRODUCT INFORMATION

REF CODE	WIDTH	LENGTH
SSN-002-S	2 in (5.0 cm)	1.5 m
SSN-002-L	2 in (5.0 cm)	10 m
SSN-003-S	3 in (7.5 cm)	1.5 m
SSN-003-L	3 in (7.5 cm)	10 m
SSN-004-S	4 in (10.0 cm)	1.5 m
SSN-004-L	4 in (10.0 cm)	10 m
SSN-006-S	6 in (15.0 cm)	1.5 m
SSN-006-L	6 in (15.0 cm)	10 m






FEATURES

- Comfortability – smooth application with less pressure and no pain to the patient
- Conformability & mouldability – extreme-soft fabric shapes easily to the contours of the body
- Generally no gloves are necessary for handling
- Tack-free with no layer separation
- Available pre-cut and in rolls across a range of widths

APPLICATION INSTRUCTIONS

01

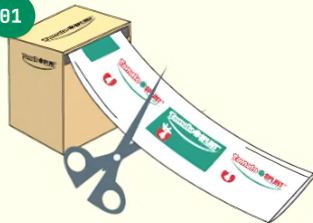

Cut to length

STEP / 01

Take Tomato Splint with its sealed foil pouch from the box, measure the length, and cut it off with bandage scissors.

02



Reseal the pack

STEP / 02

Seal the box immediately with the clip provided inside to prevent moisture from entering the package.

03




Choose the size

STEP / 03

Choose the desired size for the affected body part. Generally, no gloves are necessary.

04

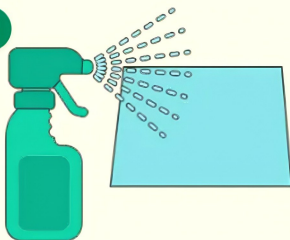


Spray

STEP / 04

Instead of soaking, spray room-temperature water on the inside (white side) of the splint for 3–5 seconds, then wipe the water off the surface.

05

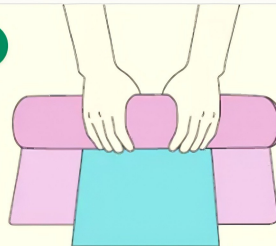


Dry off

STEP / 05

Remove excess water by laying the splint on a dry towel and rolling it up with the towel to strengthen adhesion between layers.

06

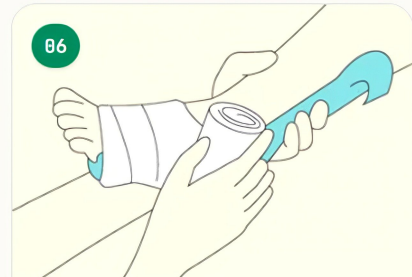

Wrap & mould

STEP / 06

Wrap Tomato Splint with a roll of elastic bandage to secure it and mould it freely to the limb.

PRODUCT INFORMATION

REF CODE	WIDTH	LENGTH	TYPE
TPS-210	2 in (5.0 cm)	10 in (25.0 cm)	Pre-cut
TPS-312	3 in (7.5 cm)	12 in (30.0 cm)	Pre-cut
TPS-335	3 in (7.5 cm)	35 in (87.5 cm)	Pre-cut
TPS-415	4 in (10.0 cm)	15 in (37.5 cm)	Pre-cut
TPS-430	4 in (10.0 cm)	30 in (75.0 cm)	Pre-cut
TPS-530	5 in (12.5 cm)	30 in (75.0 cm)	Pre-cut
TPS-545	5 in (12.5 cm)	45 in (112.5 cm)	Pre-cut
TPS-630	6 in (15.0 cm)	30 in (75.0 cm)	Pre-cut
TPS-645	6 in (15.0 cm)	45 in (112.5 cm)	Pre-cut
TRS-002	2 in (5.0 cm)	5 yds (4.5 m)	Roll
TRS-003	3 in (7.5 cm)	5 yds (4.5 m)	Roll
TRS-004	4 in (10.0 cm)	5 yds (4.5 m)	Roll
TRS-005	5 in (12.5 cm)	5 yds (4.5 m)	Roll
TRS-006	6 in (15.0 cm)	5 yds (4.5 m)	Roll
TRS-008	8 in (20.0 cm)	5 yds (4.5 m)	Roll





A VENTURE BY



BLACKCHIP IMPEX PVT. LTD., THE PARENT COMPANY OF TOMATO MEDICAL & CHEMICAL INDIA, WAS FORMED TO ACT AS A BRIDGE BETWEEN GLOBAL INNOVATION AND LOCAL EXCELLENCE: WORLD-CLASS HEALTHCARE TECHNOLOGY, VETTED BY A CLINICIAN, DELIVERED BY A GLOBAL SOURCING SPECIALIST.



Dr. Dilip Shah

The Clinical Compass.

With 35+ years of frontline experience, Dr. Shah is the guardian of our quality. He vets every product through the eyes of a clinician, ensuring that safety, precision, patient comfort and outcomes are never compromised.



Rohan Jain




The Global Architect of Sourcing.

A forward-thinking entrepreneur with a Master's in Global Family Managed Business from SP Jain, Rohan brings a sophisticated blend of heritage and innovation to the venture. His roots lie in a successful family business legacy established by his father, Rajesh Jain, where he learned the intricacies of trade and leadership, and built the strategic foundation necessary to scale.



Blackchip Impex Pvt Ltd brings
worldclass cutting edge technology to life
.
Now available in India.

Tomatomncindia.com

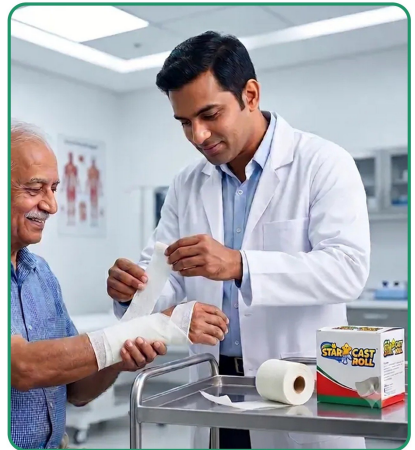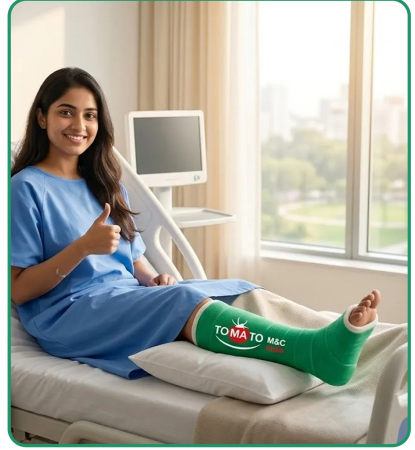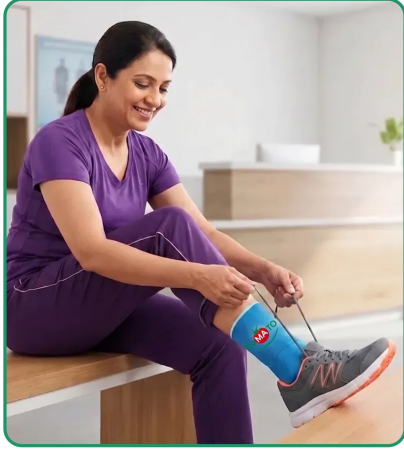

Customer care number – +91 98331 16680
Email id – tomatomncindia@gmail.com